

# MONOSWIFT®

Monofilament Poly(glycolide-co-caprolactone)  
(PGCL 25) synthetic absorbable suture

## CHARACTERISTICS

- Excellent holding strength
- Monoswift smooth tissue passage has excellent handling, exceptional suppleness, and good knot security
- Non-wicking helps prevent bacterial migration along suture line

## INDICATIONS

Monoswift® synthetic absorbable monofilament sutures are indicated for use in soft tissue approximation and/or ligation, but not for use in cardiovascular or neurological surgery, microsurgery, or ophthalmic surgery.

## COLOR

Violet or Undyed


## TENSILE STRENGTH

Implantation studies in animals indicate that PGCL suture retains approximately 70% to 80% of its original tensile strength at 7 days post implantation and approximately 40% of its original tensile strength at 14 days post implantation.

## PACKAGING

Monoswift® is available in USP sizes 6-0 through 2 (metric sizes 0.7-5).

Monoswift® is available with needles in one dozen boxes. Three dozen suture boxes available on customer request.

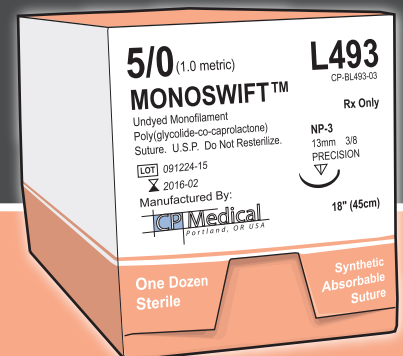

803 NE 25th Avenue, Portland, OR 97232 USA  
ph. 800 950-2763 • fax 503 230-9993 • www.cpmedical.com